

HBS ENTERPRISES

IRT-6

Infra-Red Device

Operators Manual



April 2005

HBS-IRT-6 Rev 3.0

HBS Enterprises Inc.

IRT-6
Infra-Red Unit

Operators Manual

The information in this manual is proprietary to HBS Enterprises Inc. and may not be reprinted or otherwise copied without the express written permission of HBS Enterprises Inc.

©2005 HBS Enterprises, Inc.

Table of Contents

Basic Operation of the IRT-6	4
Output Settings	5
Warranty and Return Policy.....	6
Terms of Use	7

HBS ENTERPRISES INC.
Electronic Product Design and Fabrication
3096 Hales Mill Road
Asbury, IA 52002
(563) 556-3113 phone (563) 557-8883 fax

IRT-6 is a trademark of HBS Enterprises.

Basic Operation of the IRT-6 Infra-Red Device

- 01) Connect the 12 volt DC Power adaptor into the back of the Control Unit and then plug the power adaptor into a wall outlet. (see figure 1)
- 02) Connect one or two Transducers into one or both of the transducer output jacks on the front panel. (see figure 2)
- 03) Turn on the power by moving the toggle switch on the front of the Control Unit into the upright position. (see figure 2)
- 04) The output functions are selected by pressing the Mode Button next to the LED Display. (see figure 2)
- 05) To use the audio input feature, plug the audio jacks provided into both the back of the control unit, and into a stereo, cd player, or other audio device and set the output to 0. (see figure 1) NOTE: DO NOT turn up the volume of the audio device past halfway. This could damage the IRT Device.

Figure 1


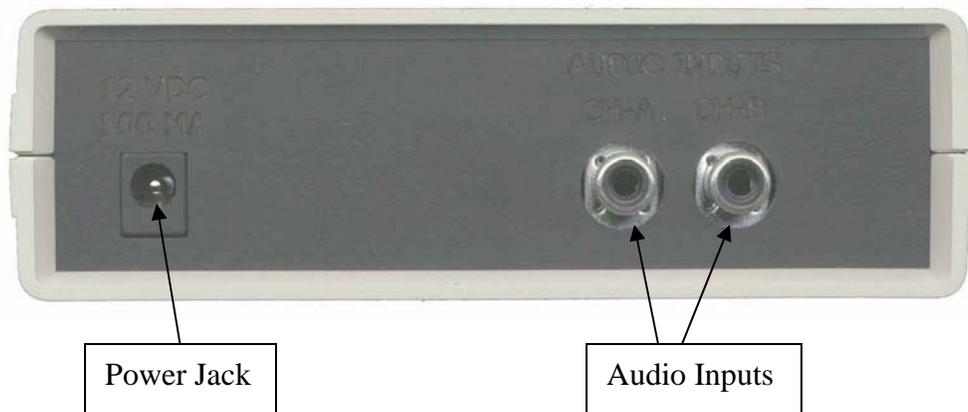


Figure 2

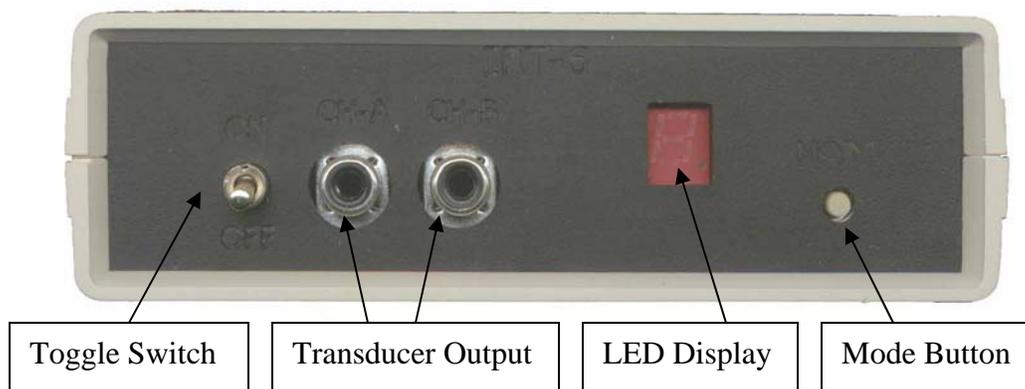

Output Settings

MODE 0: CONTINUOUS MODE / EXTERNAL AUDIO INPUT (REV-2)

MODE 1: ELF WAVE 02.28 HZ FREQUENCY SWEEPS 10-2500 HZ

MODE 2: ELF WAVE 03.42 HZ FREQUENCY SWEEPS 2500-5000 HZ

MODE 3: ELF WAVE 04.56 HZ FREQUENCY SWEEPS 72-74 HZ

MODE 4: ELF WAVE 05.70 HZ FREQUENCY SWEEPS 145-147 HZ

MODE 5: ELF WAVE 06.84 HZ FREQUENCY SWEEPS 292-294 HZ

MODE 6: ELF WAVE 07.98 HZ FREQUENCY SWEEPS 585-589 HZ

MODE 7: ELF WAVE 09.12 HZ FREQUENCY SWEEPS 1172-1176 HZ

MODE 8: ELF WAVE 10.26 HZ FREQUENCY SWEEPS 2347-2351 HZ

MODE 9: ELF WAVE 11.40 HZ FREQUENCY SWEEPS 4696-4700 HZ

Warranty and Return Policy

The IRT-6 Infra-Red Device is warranted for a period of one year after purchase against any defects or failures due to the manufacturing process. If the IRT-6 fails to operate, contact your distributor to arrange for warranty replacement or repair. If you are not satisfied with the IRT-6 within 1 week after receiving the unit, your purchase will be refunded or credit issued, at our discretion.

Disclaimer

This Document and all contents of this document are provided on an “as is” basis without warranties of any kind, either expressed or implied, including without limitation warranties of title or implied warranties of merchantability or fitness for a particular purpose. You acknowledge, by your use of the IRT-6, that your use of the above said product is at your sole risk. That you assume full responsibility for all costs associated with returning the IRT-6 for warranty and other repairs or servicing, and that HBS Enterprises Inc or its officers shall not be responsible or liable for any damages of any kind related to your use of this product.

Inaccuracy Disclaimer

From time to time there may be information in this document that contains typographical errors, inaccuracies, or omissions that may relate to product descriptions, pricing, and availability. We reserve the right to correct any errors, inaccuracies or omissions and to change or update information at any time without prior notice (including after you have submitted your order). We apologize for any inconvenience this may cause you. If you are not completely satisfied with your IRT-6, please see our return policy as described above.

Indemnification

You agree to defend, indemnify and hold HBS Enterprises Inc and officers harmless from and against any and all claims, damages, costs and expenses, including attorneys’ fees, arising from or related to your use of the IRT-6.

Miscellaneous

Unless otherwise specified and except to the extent that HBS Enterprises Inc products are offered for sale worldwide through our website, our website and the contents thereof are displayed solely for the purpose of promoting HBS Enterprises Inc products available in Canada and the United States and select foreign markets.

Terms of Use

HBS ENTERPRISES INC. does not suggest, Endorse, or imply the fact that in anyway that HBSIRT-6 is a treatment or cure for any ailment or disease. We simply sell this device for research and experimentation by those persons researching the effects 850nm pulsed infrared LED's at different frequencies. This information in no way should replace your personal physician's advice. The information contained herein is not medical advice and is not intended to replace the advice or attention of health care professionals. Nor is this information intended to act as a "prescription" for treatment. We advise you to consult your health care provider before proceeding with any type of self directed research program.

HBS ENTERPRISES INC. makes no representations concerning the efficacy, appropriateness, or suitability of any of its products or treatments. Use this device at your own risk. In view of the possibility of human error, neither the author, HBS ENTERPRISES INC nor any other party involved in providing this document, warrant that the information contained therein is in every respect accurate or complete. We are not responsible, nor liable for any errors or omissions that may be found in this document or for the results obtained from the use of such information. The entire risk as to use of this document and the device are assumed by the user. This information is for educational purposes only. If you are ill, see a health care professional. These statements have not been evaluated by the Food and Drug Administration.

However, it is your God-given right and your constitutional right under the right of privacy of the Ninth Amendment of the United States Constitution (See *Griswold vs. Connecticut* 381 US 479, June 7, 1965) to prescribe treatment for yourself, however this can involve risk. If you choose to use the device and information in this document without the approval of a health professional, you must assume the risk.